

# Zirkonzahn®

*Human Zirconium Technology*



**ALL ZIRCONIA**

## *Restoration on six implants*


Initial situation: The treatment plan for this patient with severe overbite was to move the entire arch towards dorsal. In case of extreme retro-positioning like this, the use of a coloured zirconia frame was indicated as it simplifies the reconstruction with screw access channels located buccally in front of the teeth. The latter was veneered with pink *ICE-Zirconia-Ceramic-Tissue* material. The screw access channels were plugged with tissue-coloured composite (Gradia, GC).

Unfortunately, two of the screw channels ended up in the canines' area and both had to be capped by cement-on crowns.



## *Restoration on six implants*

Dental Laboratory of Antonio Lazetera, Italy - Dr. Andrea Damiola, Italy

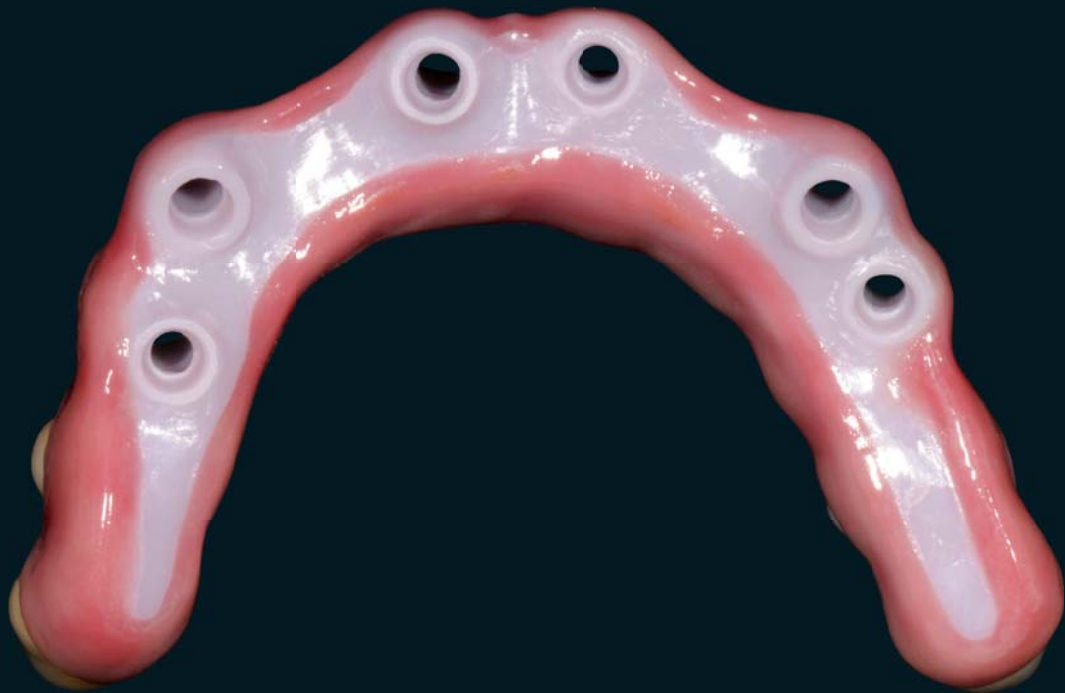


# Zirkonzahn®



## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)

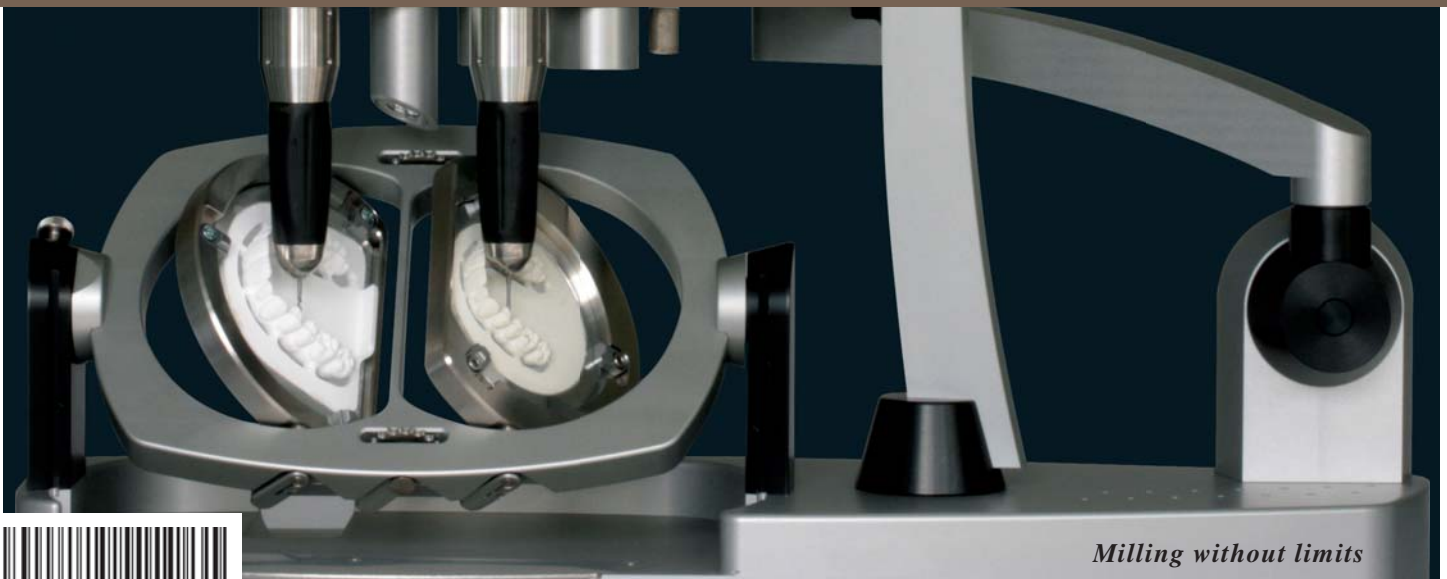

WEAA5411=

*Milling without limits*