

# Zirkonzahn®

*Human Zirconium Technology*



**100% PRETTAU ZIRCONIA**

*“Perfection through perseverance” Enrico Steger*



*Zirkonzahn World Wide*



*Zirkonzahn Germany*



*Zirkonzahn USA*



*Zirkonzahn Mexico*

## *Prettau Zirconia*

Plan for the following case was the construction of a full-anatomical-contour Prettau Bridge on six Straumann implants. Abutment angles in the areas of 17 and 27 called for primary custom abutments plus separate secondary bridge segments with lateral screws. These secondary segments were parallel-milled ('zero' degrees) and fastened to the primary bridge frame with help of the Zirkonzahn Screw-Tec-System. The four anteriors were only labially layered with veneer porcelain.

When veneering the tissue coloured areas of the restoration, special attention was paid to distinguish between gingival and alveolar crests by using varying shades of pink in order to copy the natural appearance of soft and hard tissues.



## *The Prettau Bridge*

Juergen Auffahrt – Austria

In cooperation with Master Technician Christian Hannker, Dental Laboratory Bellmann & Hannker (Germany)




# Zirkonzahn®



## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)

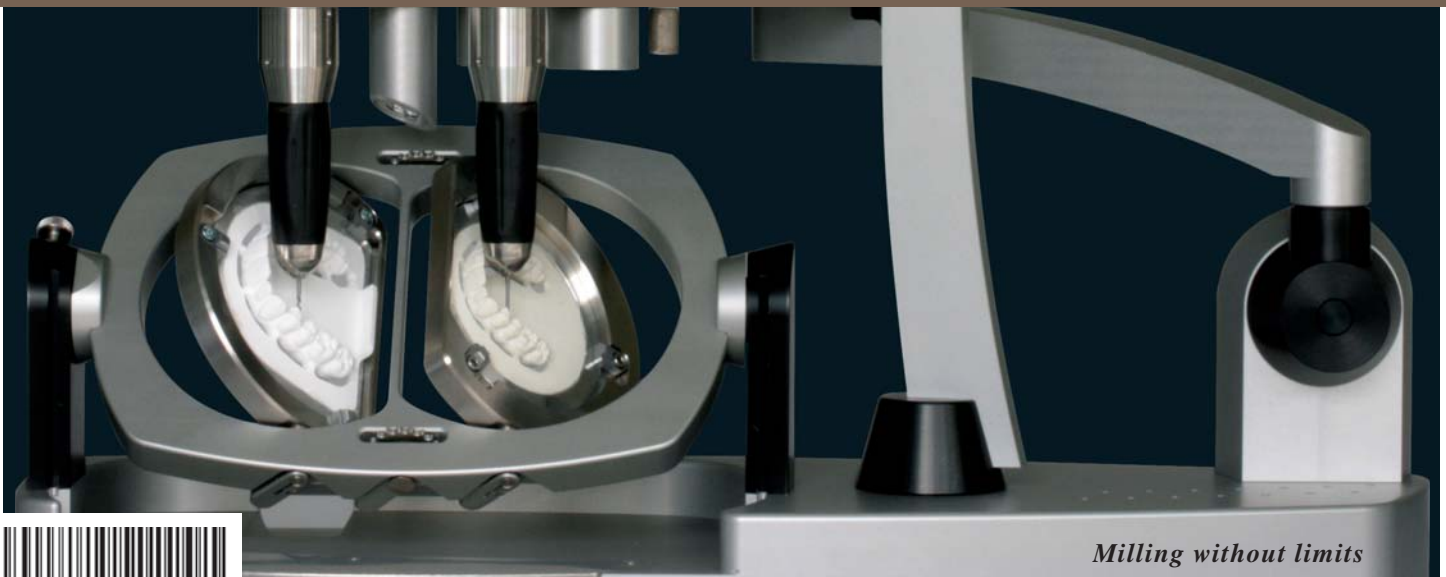

WEAA5422=

*Milling without limits*