

Zirkonzahn[®]

Human Zirconium Technology




100% PRETTAU[®] ZIRCONIA

A lower patient-removable full denture – guaranteed to last

The idea was to produce a bar-born removable lower restoration made entirely from ceramic material. Two zirconia bars were milled; both fitted with slide-attachments, one at each end. Through Teflon-sleeve matrixes over the attachments sufficient friction is achieved to hold the denture in place and give the desired stability. Such longevity is guaranteed.



The bars are milled by reducing a full-anatomical-contour mock-up frame; in this way optimal bar position within the restoration is ensured. The entire prosthesis is made from 100% full zirconia. Only the pink tissue was layered with veneer porcelain. The result is a Prettau® Bridge, worn over bars and entirely patient-removable.

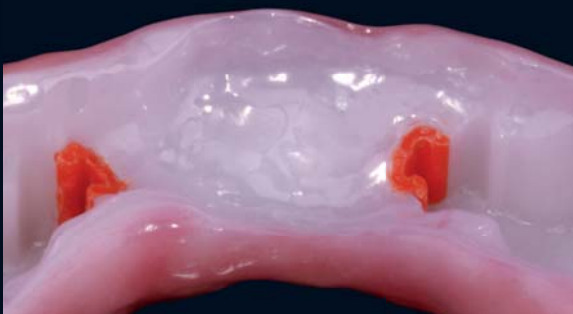

The Prettau® Bridge

Georg Walcher, Zirkonzahn Italy

Dr. Nico Bühler

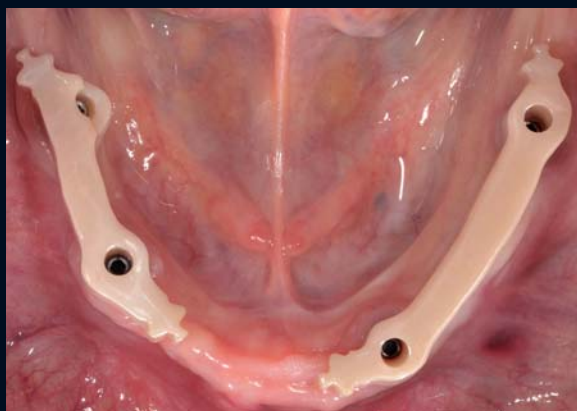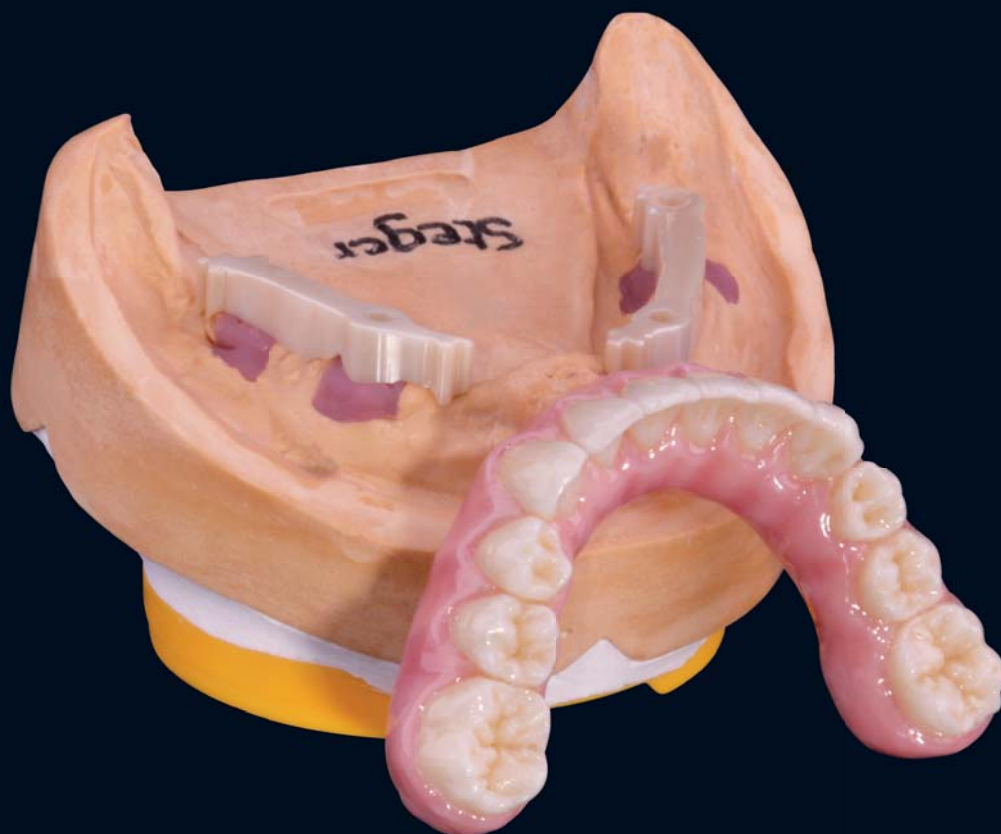


Zirkonzahn®



HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN Worldwide – Tel +39 0474 066 680 – info@zirkonzahn.com – www.zirkonzahn.com

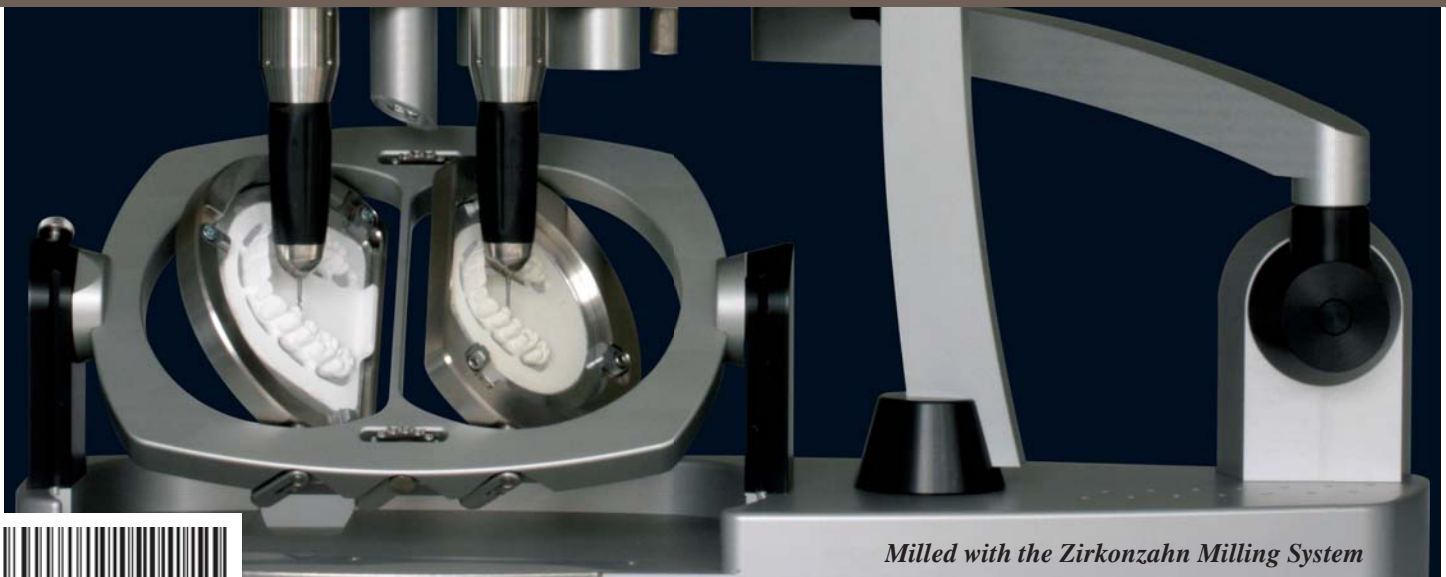

Milled with the Zirkonzahn Milling System



WEAA5209=